

C-track

Overhead tracking: Go to the wall

Primary application

Home care
Nursing homes
Hospital

Ergolet's C-track is the newest member of the Luna tracking family, ready to carry on the proud tradition of innovative Scandinavian design.

This unique wall mounted track solution is in reality a thin, flat panel offering a discreet and unobtrusive alternative to unsightly wall or ceiling brackets that are conventionally used. Apart from being an attractive alternative, C-track is simple and quick to instal making it extremely cost effective, which is why C-track is fast becoming a preferred choice for home and institutional settings.

Flexibility is a key feature of C-track, allowing traverse systems to be transferred simply, easily and quickly from one room to another. Solve the question of deciding which rooms need overhead tracking: C-track is the answer.

C-track uses the same channel to accommodate both the moving carts and the wall fixations and, when not in use, this channel may be concealed using the Ergolet cover plate transforming C-track into a fully covered, discreet wall panel.

Additionally, C-track has a narrow channel on the underside, a practical feature for hanging light book shelves, pictures or other wall decorations without having to drill holes or knock nails into walls.

Key standard functions

Can be mounted on almost any type of wall construction

Accessories


A range of accessories are available for C-track such as brackets for hanging pictures etc.

A range of suitable brackets for different installation purposes are also available.



PRODUCT HIGHLIGHTS


- Class leading design and functionality.
- Can be mounted on almost any type of wall construction.
- Quick and simple to install. Cost saving.
- Telescopic cart to compensate for varying wall widths.
- Profile can be used to support light room dividers or curtains.
- Clean and simple state-of-the-art-design.

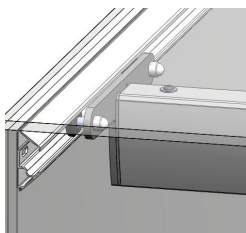


Unique telescopic cart to compensate for varying wall widths.




Quick and easy install, Saves time and reduces costs.

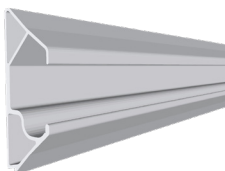

Discreet thin profile.

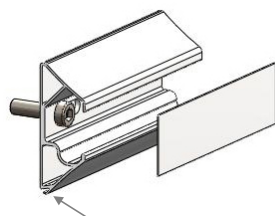

The profile has an extra track underneath, which can be used to hang pictures, book shelves or other wall decorations thus eliminating unsightly holes in the wall.


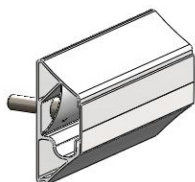


If not in use, C-track easily becomes part of the room decoration when a cover plate is clicked into place.

Version	Item-number
C-track 4,3m	016-00990
C-track 5,5m	016-00991
C-track 7m	016-00992

TECHNICAL SPECIFICATION

Load on profile	200kg and 275kg
Measurements of profile	33 mm x 87 mm
Weight	1,14 kg/m
Material	Aluminium
Mounting examples	Plasterboard, wood, concrete, brick
Colour	White powder coating

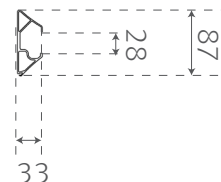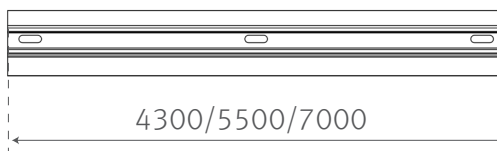

Consult the Technical manual for further specifications. All figures and measurements are for guidance only and may differ in actual situations. Ergolet reserves the right to changes without notice.

Product photos may include accessories. Consult your dealer.



date 06.01.2012 - ver. 5

WinnCare Nordic ApS
 Taarnborgvej 12C
 4220 Korsoer
 Denmark
 tel. +45 70 27 37 20
 e-mail: info@winncare.dk

www.winncare.dk

WINNCARE
 N o r d i c